

## Key Features

### DEHP / Plasticiser free

No plasticiser containing plastic anywhere in fluid / drug path – no risk of plasticiser induced toxicity

### Latex free

Reduces the chance of a severe allergic reaction

### Negligible common fluid space

Patented multilumen hub keeps infusions separate until distal tip of connector so helping to prevent drug / drug reactions

### Lipid resistant plastics

Safe for use with lipid containing drugs e.g. propofol

### Fully rotating male Luer collar

'Self locates' on locking and prevents line spaghetti

## Gravity Line

### Kink resistant gravity line

Nylon braiding built into catheter wall helps maintain free-flowing drip

### Low sorbing plastic

Polyurethane catheter prevents changes in bioavailability of infused drug

### Low resistance back check valve – opens with minimal pressure

Effectively prevents back flow of drugs or blood into gravity line

## Drug Line(s)

### Low sorbing plastic

Trilayer co-extruded catheter with polyethylene inner lining prevents changes in bioavailability of infused drug

### Normally dosed anti-siphon valve(s)

One way valve opens with higher crack pressure, preventing accidental overdose from siphoning from pump or syringe

### Colour coded pinch clamp

Prevents accidental bolus delivery during syringe handling

### Detachable long line extensions on drug line(s)

Permits exchange of pumped infusion for new drug e.g. at end of operation

## A2-Set


Ref: 7200A

### 1 gravity ('drip') line

- braided
- female Luer back check valve
- Ø: 3.70 x 1.85 mm
- L: 8 cm
- Vol: 0.5 ml

### 1 infusion line

- white female Luer anti-siphon valve
- white pinch clamp
- Ø: 3.10 x 1.70 mm
- L: 10 cm
- Vol: 0.6 ml

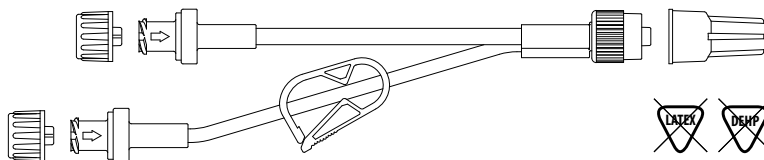

## A2L-Set


Ref: 7200L

### 1 gravity ('drip') line

- braided
- female Luer back check valve
- Ø: 3.70 x 1.85 mm
- L: 8 cm
- Vol: 0.5 ml

### 1 infusion line

- white female Luer anti-siphon valve
- white pinch clamp
- detachable extension
- Ø: 3.10 x 1.70 mm
- L: 210 cm
- Vol: 5.1 ml

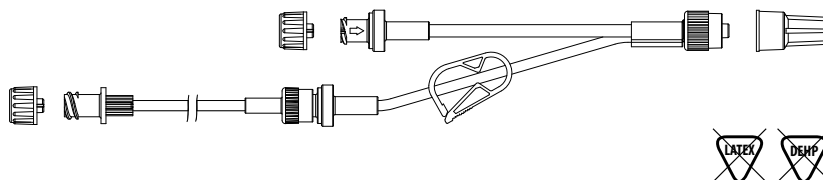

## A22BL-Set


Ref: 7220BL

### 1 gravity ('drip') line

- braided
- female Luer back check valve
- Ø: 3.70 x 1.85 mm
- L: 8 cm
- Vol: 0.5 ml

### 1 infusion line

- white female Luer anti-siphon valve
- white pinch clamp
- bonded (non-detachable) extension
- Ø: 3.10 x 1.70 mm
- L: 200 cm
- Vol: 4.75 ml

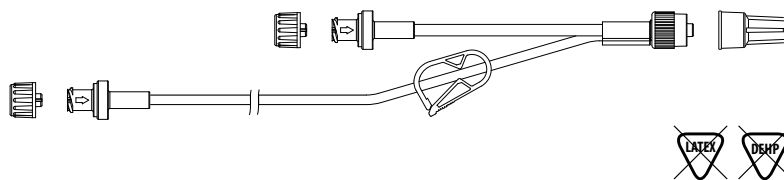

## A2NFV-Set


Ref: 7200ANF

### 1 gravity ('drip') line

- braided
- female Luer back check valve
- Ø: 3.70 x 1.85 mm
- L: 8 cm
- Vol: 0.5 ml

### 1 infusion line

- needle-free valve
- white pinch clamp
- Ø: 3.10 x 1.70 mm
- L: 10 cm
- Vol: 0.6 ml

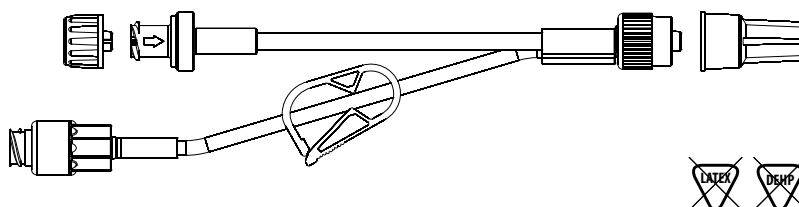

## A3-Set

Ref: 7300A

### 1 gravity ('drip') line

- braided
  - female Luer back check valve
- Ø: 3.70 x 1.85 mm  
L: 8 cm  
Vol: 0.5 ml


Ø: 3.10 x 1.70 mm  
L: 10 cm  
Vol: 0.6 ml

### 1 infusion line

- white female Luer anti-siphon valve
- white pinch clamp
- bonded (non-detachable) extension

### 1 infusion line

- green female Luer anti-siphon valve
  - green pinch clamp
  - bonded (non-detachable) extension
- Ø: 3.10 x 1.70 mm  
L: 12 cm  
Vol: 0.3 ml

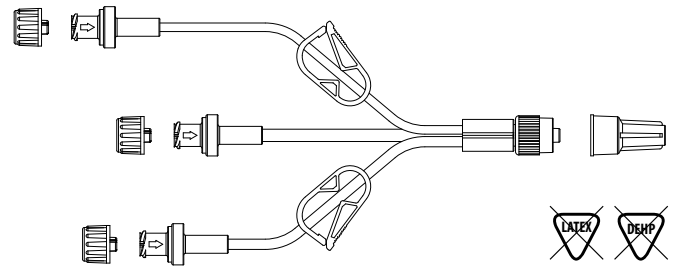

## A3L-Set

Ref: 7300AL

### 1 gravity ('drip') line

- braided
  - female Luer back check valve
- Ø: 3.70 x 1.85 mm  
L: 8 cm  
Vol: 0.5 ml


Ø: 3.10 x 1.70 mm  
L: 185 cm  
Vol: 5.0 ml

### 1 infusion line

- white female Luer anti-siphon valve
- white pinch clamp
- detachable extension

### 1 infusion line

- green female Luer anti-siphon valve
  - green pinch clamp
  - detachable extension
- Ø: 3.10 x 1.70 mm  
L: 187 cm  
Vol: 4.7 ml

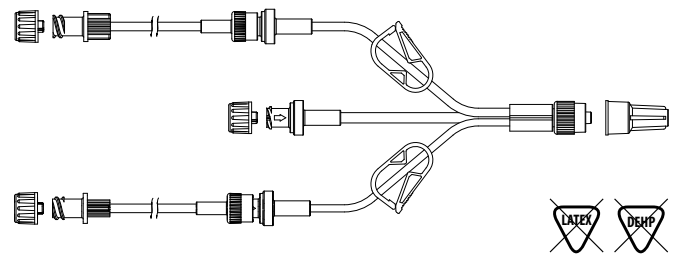

## A32BL-Set

Ref: 7320BL

### 1 gravity ('drip') line

- braided
  - female Luer back check valve
- Ø: 3.70 x 1.85 mm  
L: 8 cm  
Vol: 0.5 ml


Ø: 3.10 x 1.70 mm  
L: 200 cm  
Vol: 4.75 ml

### 1 'Siamese' infusion line

- white female Luer anti-siphon valve
- white pinch clamp
- fixed extension

### 1 'Siamese' infusion line

- green female Luer anti-siphon valve
  - green pinch clamp
  - fixed extension
- Ø: 3.10 x 1.70 mm  
L: 200 cm  
Vol: 1.5 ml

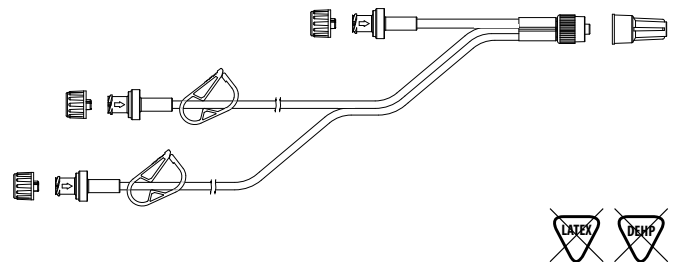

## A33BL-Set

Ref: 7330BL

### 1 gravity ('drip') line

- braided
  - female Luer back check valve
- Ø: 3.70 x 1.85 mm  
L: 8 cm  
Vol: 0.5 ml


Ø: 3.10 x 1.70 mm  
L: 300 cm  
Vol: 7.0 ml

### 1 'Siamese' infusion line

- white female Luer anti-siphon valve
- white pinch clamp
- fixed extension

### 1 'Siamese' infusion line

- green female Luer anti-siphon valve
  - green pinch clamp
  - fixed extension
- Ø: 3.10 x 1.70 mm  
L: 300 cm  
Vol: 2.5 ml

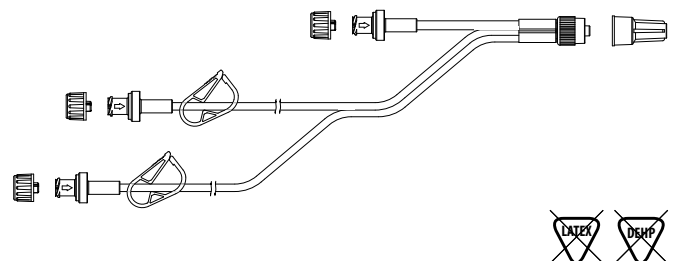

## A3NFV-Set

Ref: 7300ANF

### 1 gravity ('drip') line

- braided
  - back check valve
  - female Luer back check valve
- Ø: 3.70 x 1.85 mm  
L: 8 cm  
Vol: 0.5 ml

Ø: 3.10 x 1.70 mm  
L: 10 cm  
Vol: 0.6 ml

### 1 infusion line

- white anti-siphon valve
- female Luer needle-free valve
- white pinch clamp

### 1 infusion line

- green anti-siphon valve
  - female Luer needle-free valve
  - green pinch clamp
- Ø: 3.10 x 1.70 mm  
L: 12 cm  
Vol: 0.3 ml

