

# Key Features

## Promote infection control

No risk of bacterial colonisation in stop cock (three-way tap)

## Latex free

Reduces the chance of a severe allergic reaction

## DEHP / Plasticiser Free

No plasticiser containing plastic in fluid / drug path – no risk of plasticiser induced toxicity

## Negligible common fluid space

Patented multilumen hub keeps infusions separate until distal tip of connector so helping to prevent drug / drug reactions

## Lipid resistant plastics

Safe for use with lipid containing drugs

## Fully rotating male Luer collar

'Self locates' on locking and prevents line spaghetti

## Line clamps

Prevents accidental bolus delivery during syringe handling

## Normally closed valve

Opens after insertion of tip of syringe or male Luer from administration set

## Luer thread to proximal port

Permits screw locking of syringe or administration set to valve

## Two-way valve function


Allows aspiration of fluids as well as administration of drugs

## SafeSet 1 mini

Ref: 6100MNF

### 1 extension line with

- female Luer needle-free valve
- pinch clamp
- Ø: 2.6 x 1.0 mm
- L: 8 cm
- Vol: 0.1 ml

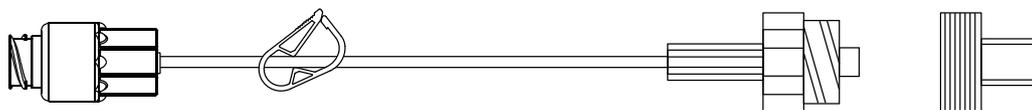


## SafeSet 2

Ref: 6200NF

### 2 extension line, each with

- female Luer needle-free valve
- pinch clamp
- Ø: 4.1 x 3.0 mm
- L: 7 cm
- Vol: 0.6 ml

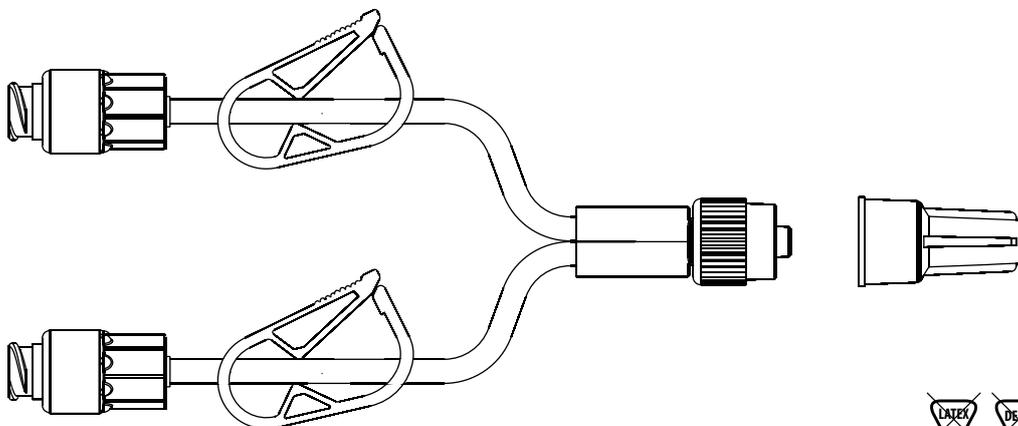


## SafeSet 2 back check

Ref: 6200BCNF

### 2 extension lines, each with

- female Luer needle-free valve
- one-way back check valve
- pinch clamp
- Ø: 4.1 x 3.0 mm
- L: 7 cm
- Vol: 0.6 ml

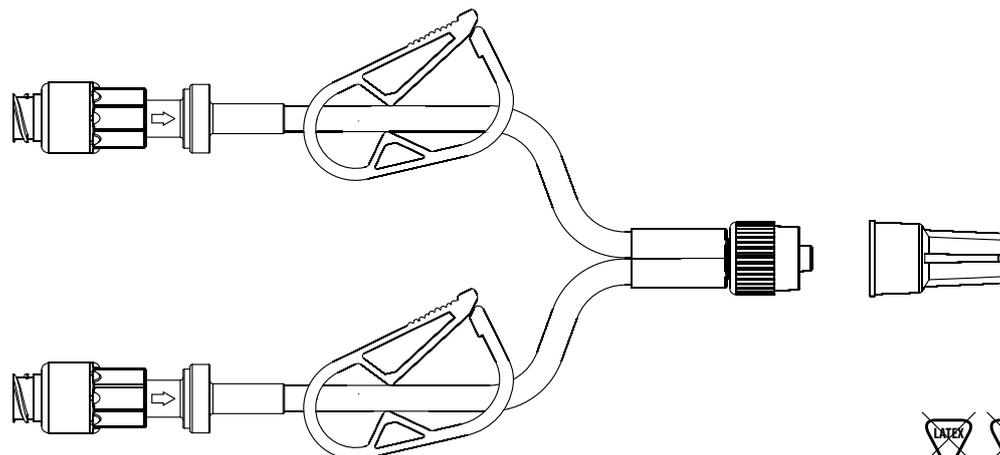


## SafeSet 2 mini

Ref: 6400MNF

2 extension polyethylene-lined mini-bore lines, each with

- female Luer needle-free valve
  - mini pinch clamp
- Ø: 2.6 x 1.0 mm  
L: 8 cm  
Vol: 0.1 ml

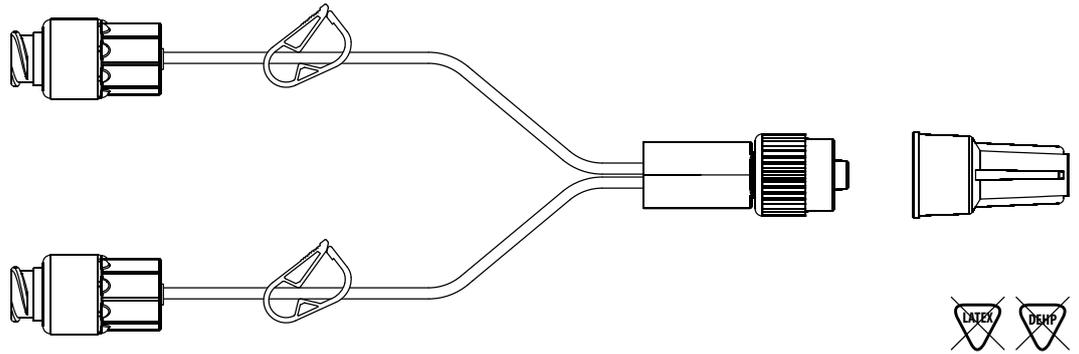


## SafeSet 3

Ref: 6300NF

3 extension lines, each with

- female Luer needle-free valve
  - pinch clamp
- Ø: 4.1 x 3.0 mm  
L: 7 cm  
Vol: 0.6 ml

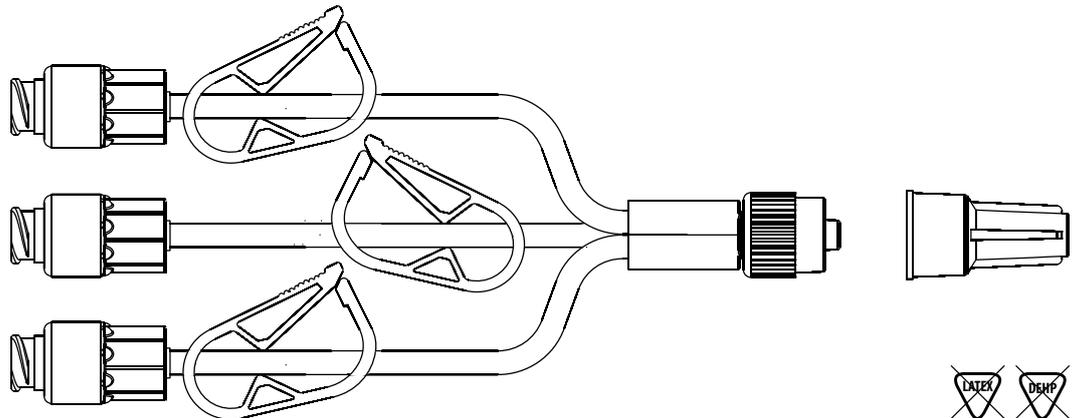


## SafeSet 3 back check

Ref: 6300BCNF

3 extension lines, each with

- female Luer needle-free valve
  - one-way back check valve
  - pinch clamp
- Ø: 4.1 x 3.0 mm  
L: 7 cm  
Vol: 0.6 ml

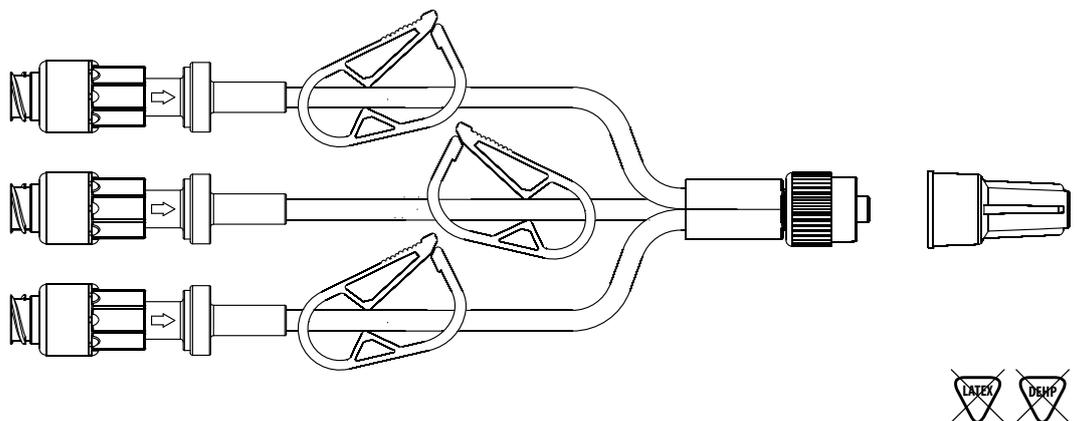

## SafeSet 3 mini

Ref: 6500MNF

3 extension polyethylene-lined mini-bore lines, each with

- female Luer needle-free valve
  - pinch clamp
- Ø: 2.6 x 1.0 mm  
L: 8 cm  
Vol: 0.1 ml

